Test your health

 **PRIMA LAB SA**
Via Antonio Monti 7
CH-6828 Balerna - SWITZERLAND
 **support@primalabsa.ch**
primalabsa.ch

19806_IFU_11_6.2 04/2025  19806-x

Registered trademarks.

INSTRUCTIONS FOR USE
GEBRAUCHSANWEISUNG
NOTICE D'UTILISATION
ISTRUZIONI PER L'USO

DRUG TEST

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PRECAUTIONS

1. The assay is designed for use with human urine only. Fresh urine specimens should be used, and the urine does not require any special handling or pretreatment.
2. Carefully read the instructions before performing the test. The test is only interpretable if the instructions are carefully respected.
3. Keep out of the reach of children.
4. The test Panel should remain in the sealed pouch until use.
5. Do not use after the expiry date printed on the label and on the protective pouch.
6. Make sure you correctly follow the "time" instructions when carrying out the test and observing the results.
7. Do not freeze. Use the test kit at temperatures between 15-30 °C.
8. After use, all components can be disposed according to local regulations; all specimen should be considered potentially hazardous and handled in the same manner as an infectious agent.
9. The test kits should be kept away from direct sunlight, moisture and heat.
10. In case of difficulties in color identification (such as Daltonism), ask for help in test reading.
11. External use only. DO NOT SWALLOW.

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
SPECIMEN HANDLING AND STORAGE

The urine specimens should be collected in a clean and dry container. Adulterants, such as bleach and/or alum, in urine specimens, may produce erroneous results regardless of the analytical method used. If adulteration is suspected, the test should be repeated with another urine specimen. Urine collected at any time of the day may be used. Urine specimens exhibiting visible precipitates should be centrifuged, filtered, or allowed to settle to obtain a clear specimen for testing. Urine specimens may be stored at 2–8 °C for up to 48 hours prior to testing. For prolonged storage, specimens may be frozen and stored below −20 °C. Frozen specimens should be thawed and mixed well before testing.

TEST PROCEDURE

Allow the test and the urine specimen, and/or controls to reach room temperature (15–30 °C) prior to testing.

- Bring the pouch to room temperature before opening it. Remove the test panel from the sealed pouch and use it within one hour.
- Remove the cap.
- With the arrow pointing toward the urine specimen, immerse the test panel vertically in the urine specimen for at least 10 to 15 seconds. Immerse the dipstick to at least the level of the wavy lines, but not above the arrow on the test panel.
- Replace the cap and place the test panel on a non-absorbent flat surface.
- Start the timer and wait for the colored line(s) to appear.

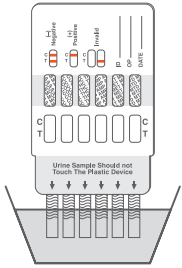RESULTS INTERPRETATION

RESULT SHOULD BE READ AT 5 MINUTES.

DO NOT INTERPRET THE RESULT AFTER 10 MINUTES.

The shade of the colored lines(s) in the Test (T) region may vary. The result should be considered negative whenever there is even a faint line. Please refer to the illustration

NEGATIVE

A colored line appears in the Control (C) region and colored lines appear in the Test (T) region. This negative result means that the concentrations in the urine sample are below the designated cut-off levels for a particular drug tested.

POSITIVE

A colored line appears in the Control (C) region and no line appears in the Test (T) region. The positive result means that the drug concentration in the urine sample is greater than the designated cut-off for a specific drug.

NON VALID

No line appears in the Control (C) region. Insufficient specimen volume or incorrect procedural techniques are the most likely reasons for Control line failure. Read the directions again and repeat the test with a new test card. If the result is still invalid, contact your manufacturer.

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Amphetamine (AMP) Amphetamine is a potent central nervous system (CNS) stimulant that is used in the treatment of attention deficit hyperactivity disorder (ADHD), narcolepsy, and obesity but it is also available on the illicit market. Amphetamines are chemically related to the human body's natural catecholamines: epinephrine and norepinephrine. Acute higher doses lead to enhanced stimulation of the central nervous system (CNS) and induce euphoria, alertness, reduced appetite, and a sense of increased energy and power. Cardiovascular responses to amphetamines include increased blood pressure and cardiac arrhythmias. More acute responses produce anxiety, paranoia, hallucinations, and psychotic behavior. The effects of amphetamines generally last 2-4 hours following use and the drug has a half-life of 4-24 hours in the body. About 30% of amphetamines are excreted in the urine in unchanged form, with the remainder as hydroxylated and deaminated derivatives. The DRUG TEST Panel yields a positive result when the concentration of amphetamines in urine exceeds the detective level of 1000 ng/mL.

Cocaine (COC) Cocaine is a potent central nervous system stimulant and a local anesthetic. Initially, it brings about extreme energy and restlessness while gradually resulting in tremors, over-sensitivity and spasms. In large amounts, cocaine causes fever, unresponsiveness, difficulty in breathing and unconsciousness.

Cocaine is often self-administered by nasal inhalation, intravenous injection and free-base smoking. It is excreted in the urine in a short time primarily as benzoylecgonine. Benzoylecgonine, a major metabolite of cocaine, has a longer biological half-life (5-8 hours) than cocaine (0.5-1.5 hours), and can generally be detected for 24-72 hours after cocaine exposure. The DRUG TEST Panel yields a positive result when the concentration of benzoylecgonine in urine exceeds the detective level of 300 ng/mL.

Marijuana (THC) THC (Δ9-tetrahydrocannabinol) is the primary active ingredient in cannabis (marijuana). When smoked or orally administered, THC produces euphoric effects. Users have impaired short-term memory and slowed learning. They may also experience transient episodes of confusion and anxiety. Long-term, relatively heavy use may be associated with behavioral disorders. The peak effect of marijuana administered by smoking occurs in 20-30 minutes and the duration is 90-120 minutes after one cigarette. Elevated levels of urinary metabolites are found within hours of exposure and remain detectable up to 4 days after smoking. The main metabolite excreted in the urine is 11-nor-Δ9-tetrahydrocannabinol-9-carboxylic acid (THC-COOH). The DRUG TEST Panel yields a positive result when the concentration of THC-COOH in urine exceeds the detective level of 50 ng/mL.

Methadone (MTD) Methadone is a narcotic analgesic prescribed for the management of moderate to severe pain and for the treatment of opiate dependence (Heroin, Vicodin, Percocet, Morphine). The pharmacology of oral methadone is very different from IV methadone. Oral methadone is partially stored in the liver for later use. IV methadone acts more like heroin. In most states you must go to a pain clinic or a methadone maintenance clinic to be prescribed methadone. Methadone is a long acting pain reliever producing effects that last from 3-5 days. Ideally, methadone frees the client from the pressures of obtaining illegal heroin, from the dangers of injection, and from the emotional roller coaster that most opiates produce. Methadone, if taken for long periods and at large doses, can lead to a very long withdrawal period. The withdrawals from methadone are more prolonged and troublesome than those provoked by heroin cessation, yet the substitution and phased removal of methadone is an acceptable method of detoxification for patients and therapists. The DRUG TEST Panel yields a positive result when the concentration of methadone in urine exceeds the detective level of 300 ng/mL.

Methamphetamine (MET) Methamphetamine is an addictive stimulant drug that strongly activates certain systems in the brain. Methamphetamine is closely related chemically to amphetamine, but the central nervous system effects of methamphetamine are greater. Methamphetamine is made in illegal laboratories and has a high potential for abuse and dependence. The drug can be taken orally, injected, or inhaled. Acute higher doses lead to enhanced stimulation of the central nervous system and induce euphoria, alertness, reduced appetite, and a sense of increased energy and power. Cardiovascular responses to methamphetamine include increased blood pressure and cardiac arrhythmias. More acute responses produce anxiety, paranoia, hallucinations, psychotic behavior, and eventually, depression and exhaustion. The effects of methamphetamine generally last 2-4 hours and the drug have a half-life of 9-24 hours in the body. Methamphetamine is excreted in the urine primarily as amphetamine, and oxidized and deaminated derivatives. However, 10-20% of methamphetamine is excreted unchanged. Thus, the presence of the parent compound in the urine indicates methamphetamine use. Methamphetamine is generally detectable in the urine for 1-4 days, depending on urine pH level. The DRUG TEST Panel yields a positive result when the methamphetamine in urine exceeds the detective level of 500 ng/mL.

Morphine (MOP) Opiate refers to any drug that is derived from the opium poppy, including the natural products, morphine and codeine, and the semi-synthetic drugs such as heroin. Opioid is more general, referring to any drug that acts on the opioid receptor. Opioid analgesics comprise a large group of substances which control pain by depressing the CNS. Large doses of morphine can produce higher tolerance levels, physiological dependency in users, and may lead to substance abuse. Morphine is excreted unmetabolized, and is also the major metabolic product of codeine and heroin. Morphine is detectable in the urine for 10-55 hours after an opiate dose. The DRUG TEST Panel yields a positive result when the concentration of morphine in urine exceeds the detective level of 300 ng/mL.

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